

Figure 1-Main page unlocked.

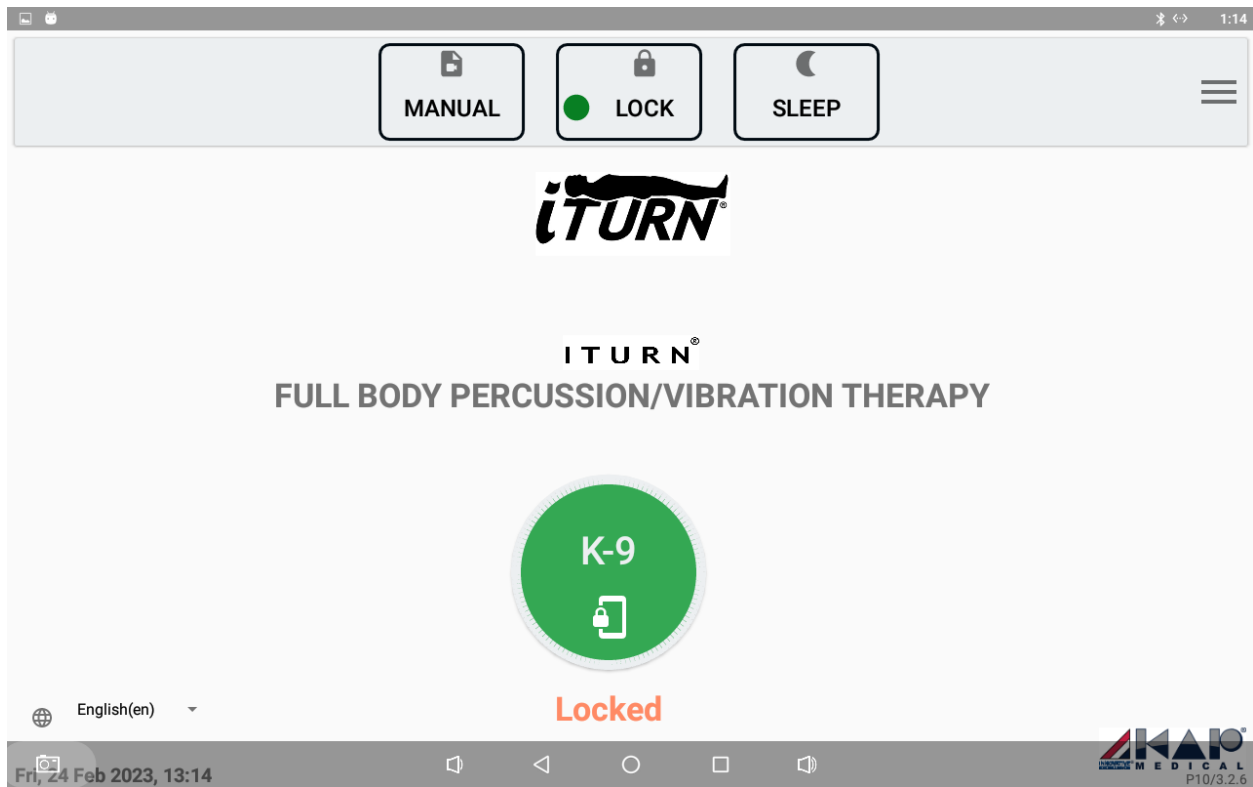


Figure 2-main page locked.

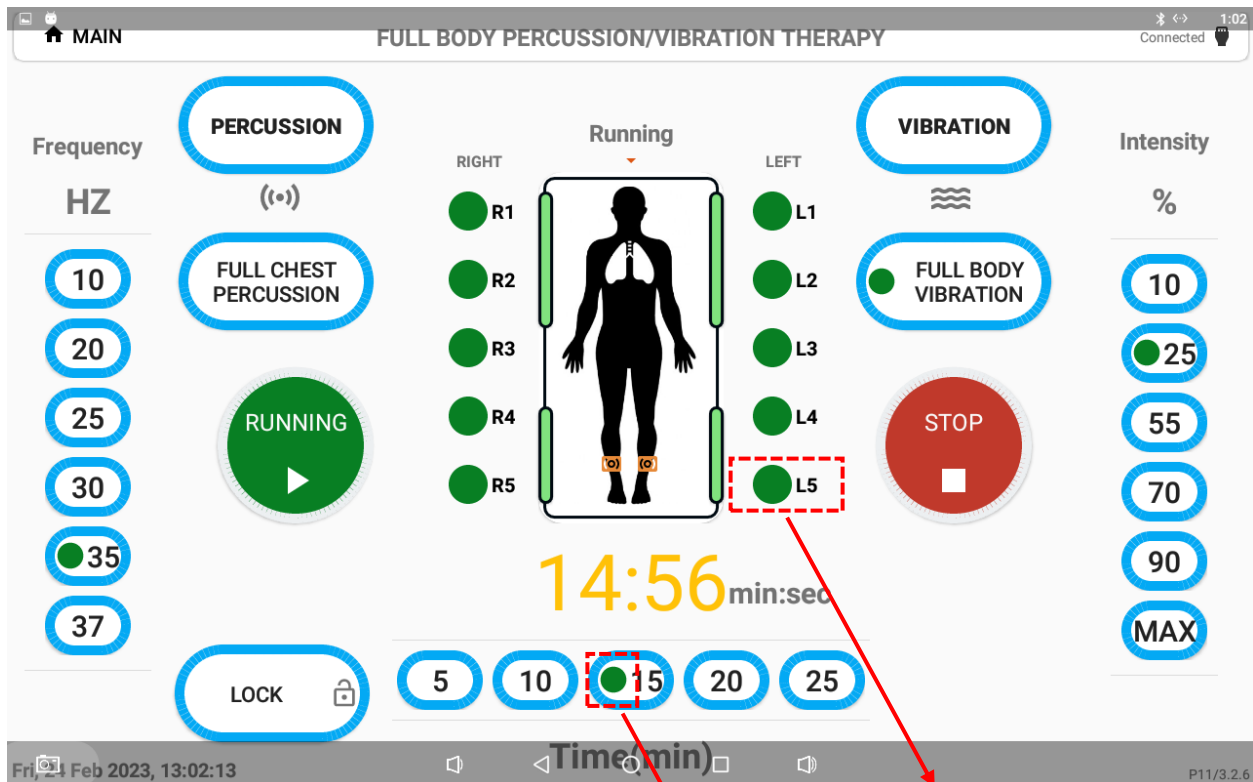


Figure 3-Therapy main page-unlocked.

Press for select or
deselect module.

Modes:

1. Percussion
2. Vibration
3. Full chest percussion
4. Full Body vibration

Parameter has been
selected when green
dot appears.

Mode-Percussion

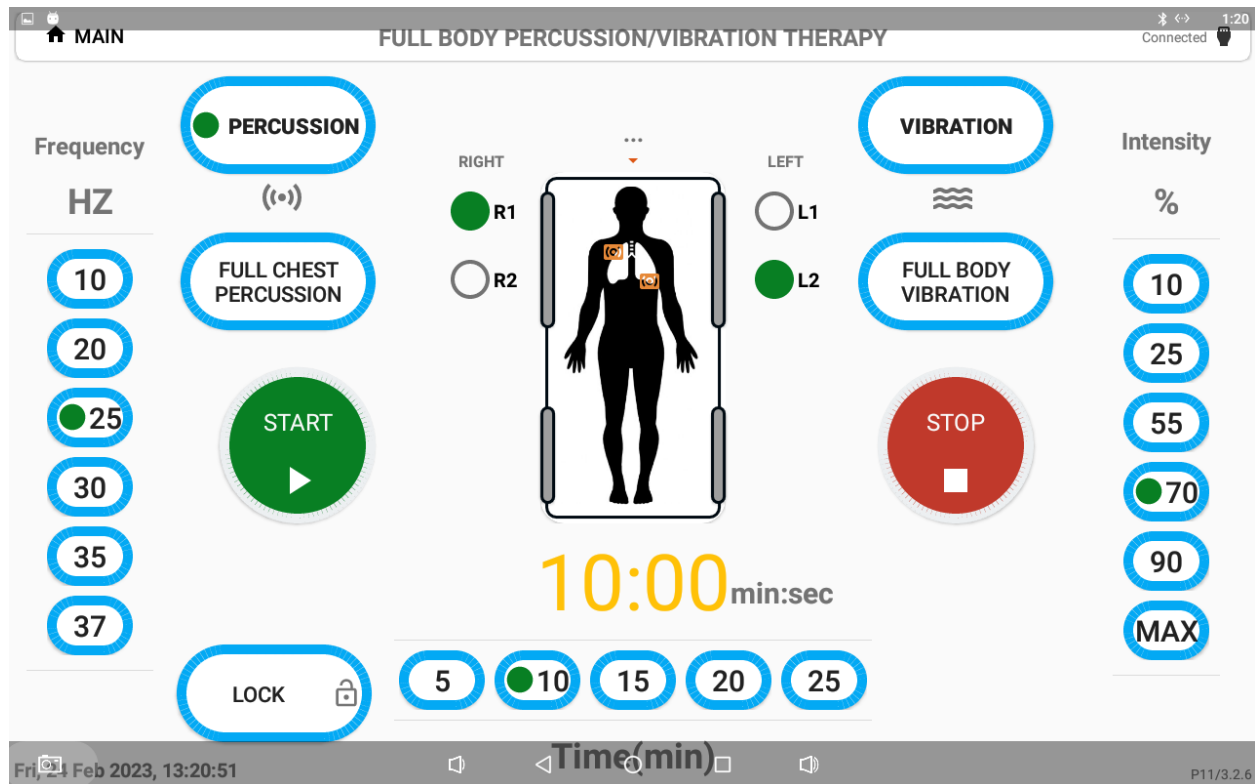


Figure 4-Mode percussion.

Select transducers options:

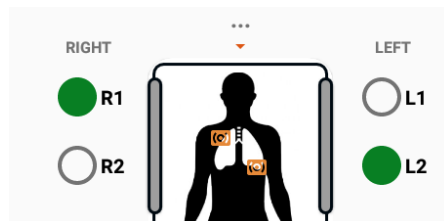


Figure 5-Mode percussion-option1

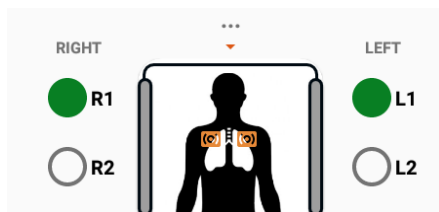


Figure 6-Mode percussion-option2

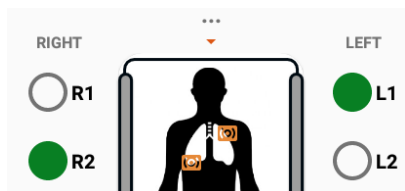


Figure 7-Mode percussion-option3.

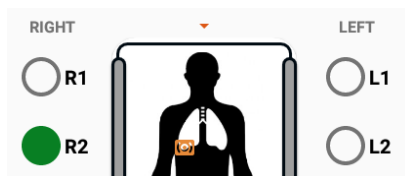


Figure 8-Mode percussion-option4.

To disable a selected module press again R# or L#

Mode – Vibration

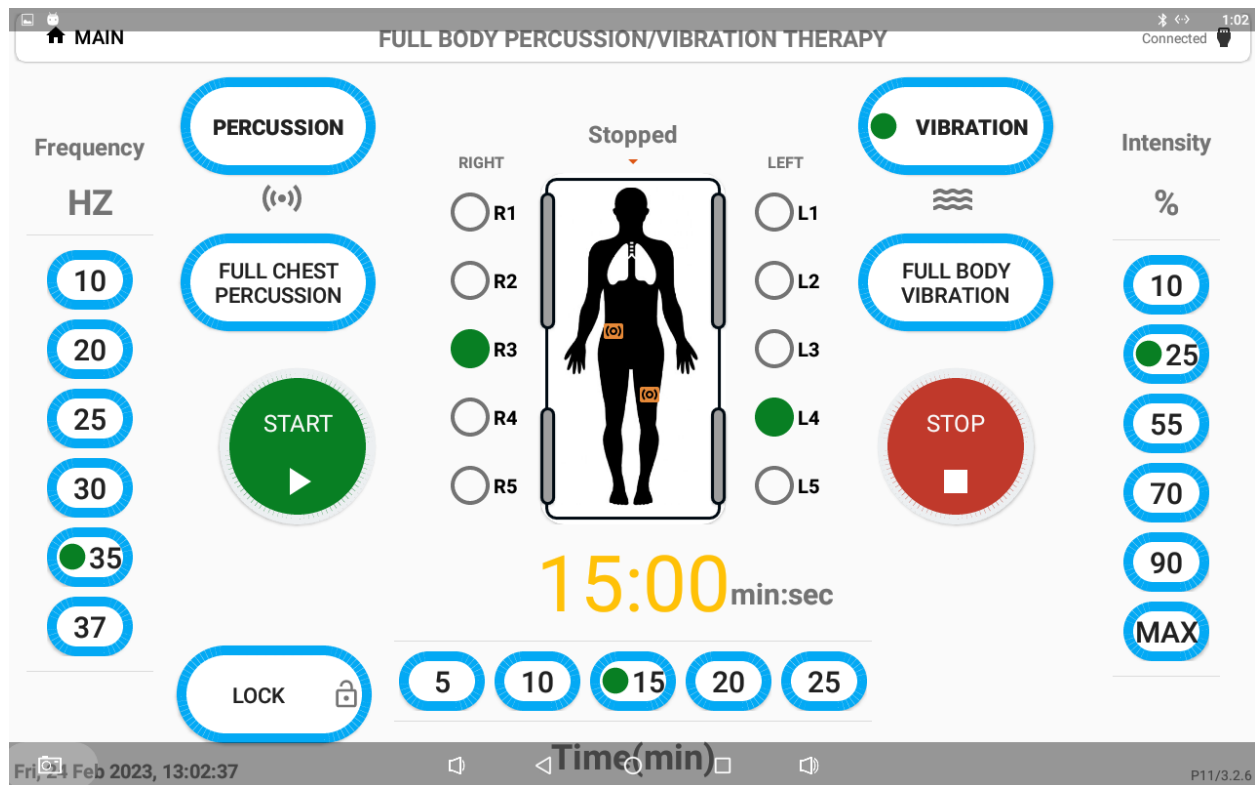


Figure 9-Mode vibration.

The modules can be selected like Percussion from 1 to 5.

Mode-Full Chest

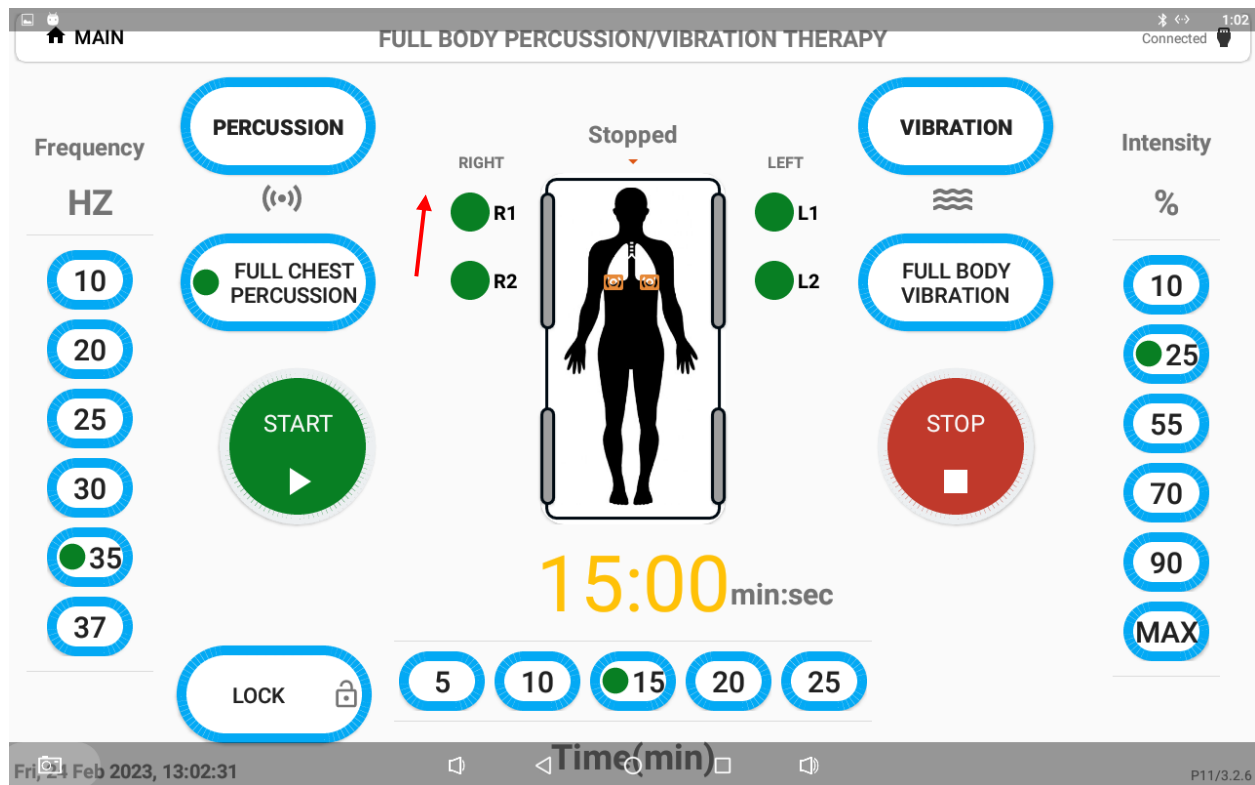


Figure 10-Mode full chest.

There are no options to select an individual module.

Mode-Full body vibration

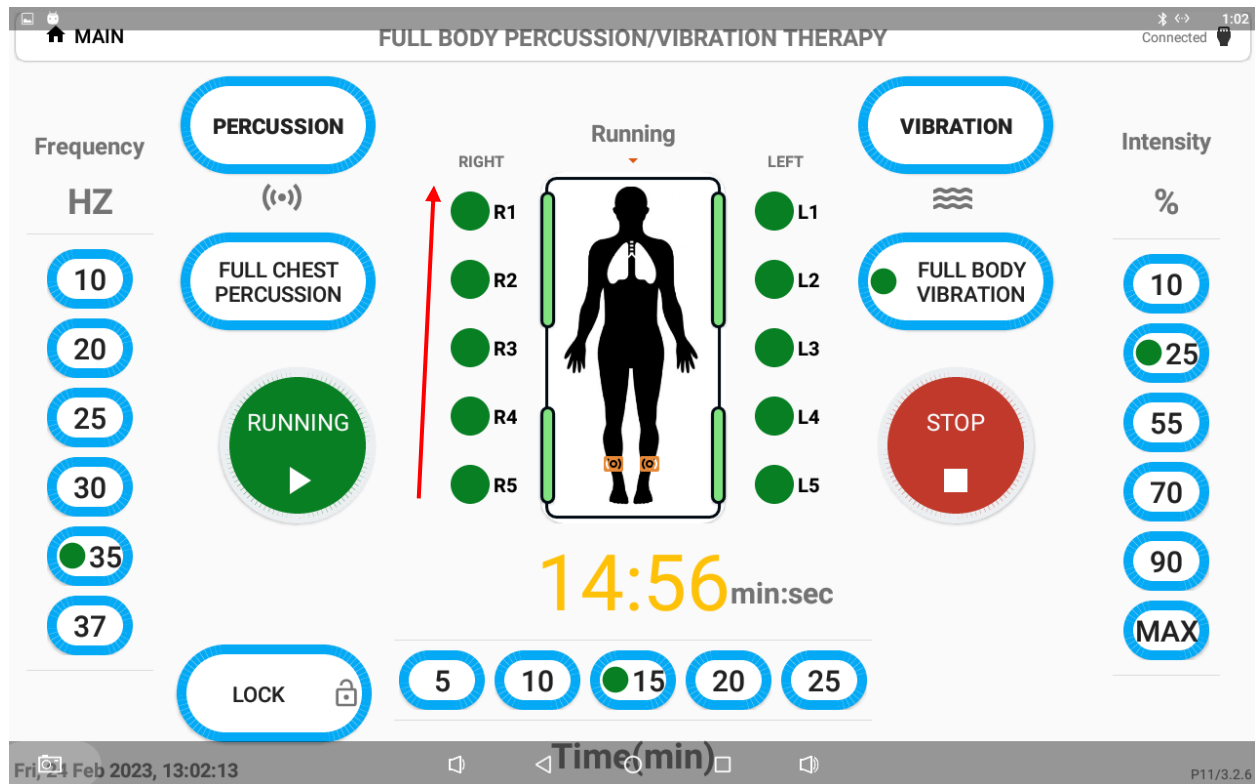


Figure 11-Mode full body vibration

The modules will run from module 5 to 1. It will repeat it until time elapsed.

Conditions before running therapy:

1. Select frequency (Figure12).
2. Select Intensity (Figure13).
3. Select time (Figure14).
4. Select mode and transducers according to mode (Figure15).
5. Side rail in lock position (Figure16).

Frequency

HZ

10

20

25

30

35

37

Figure 12-Select frequency.

Intensity

%

10

25

55

70

90

MAX

Figure 13-Select Intensity.

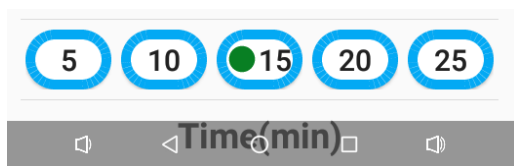


Figure 14-Select timer.

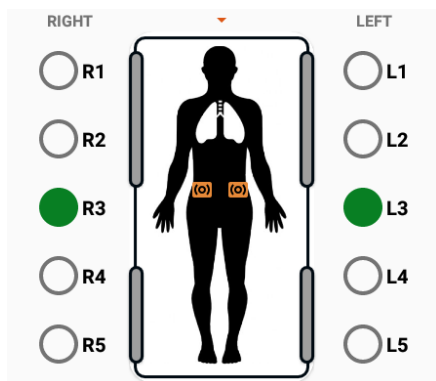


Figure 15-Select module/Transducer.


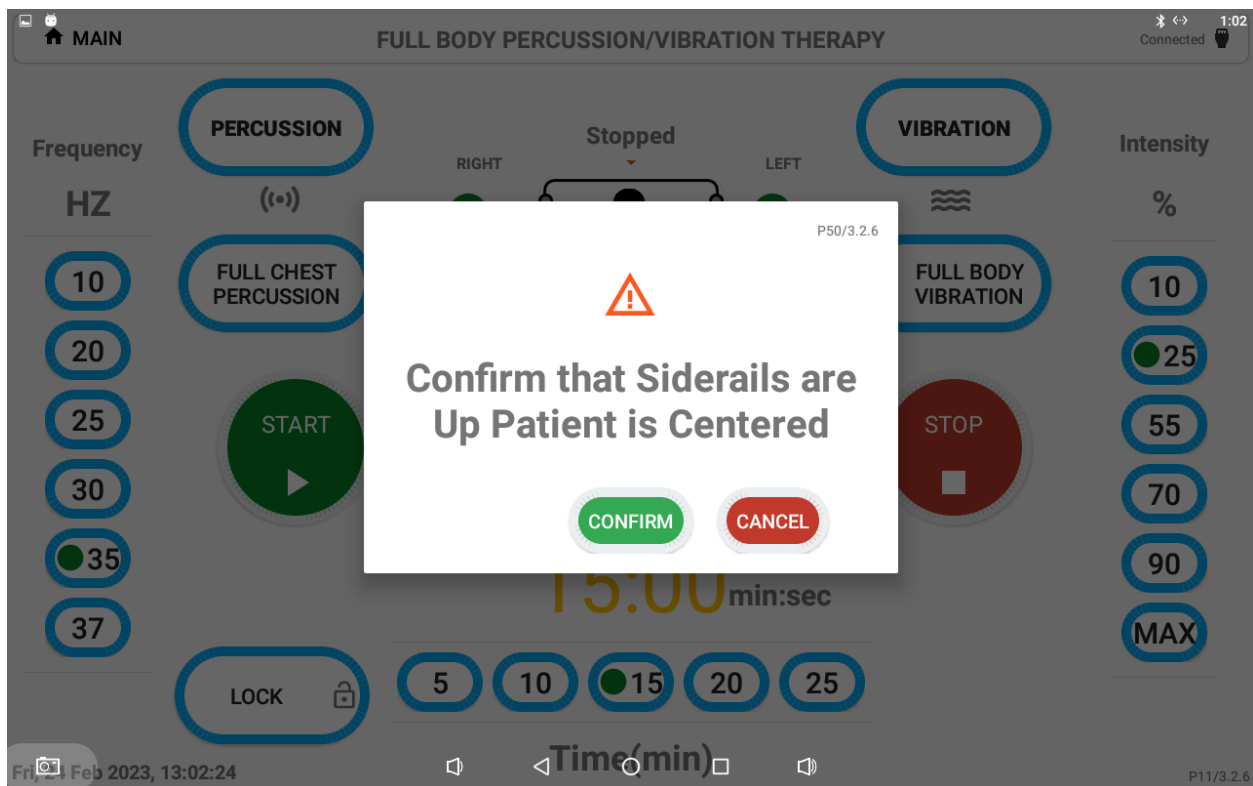


Figure 16-Confirm side rail.

Therapy running

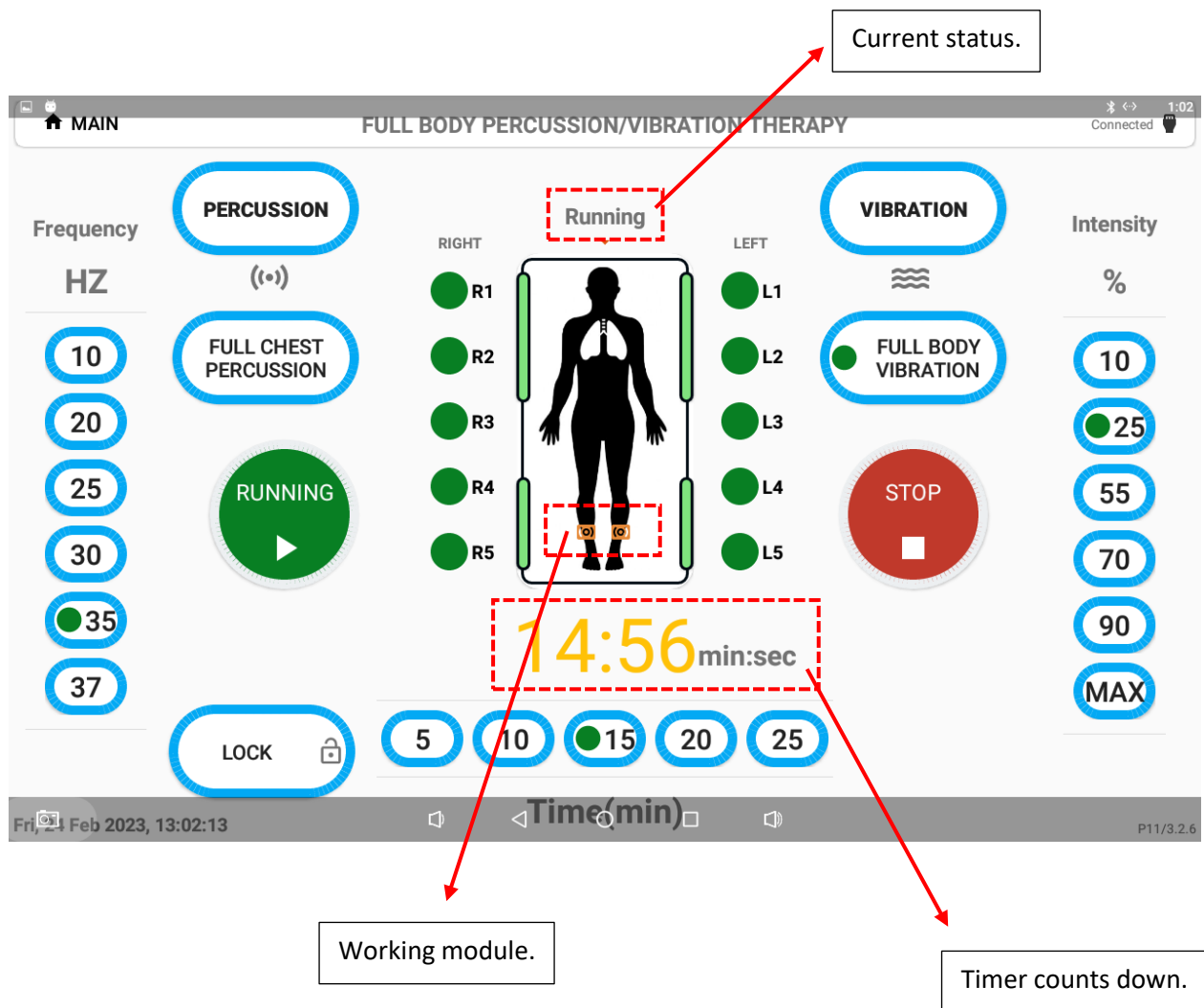

Therapy locked by nurse.

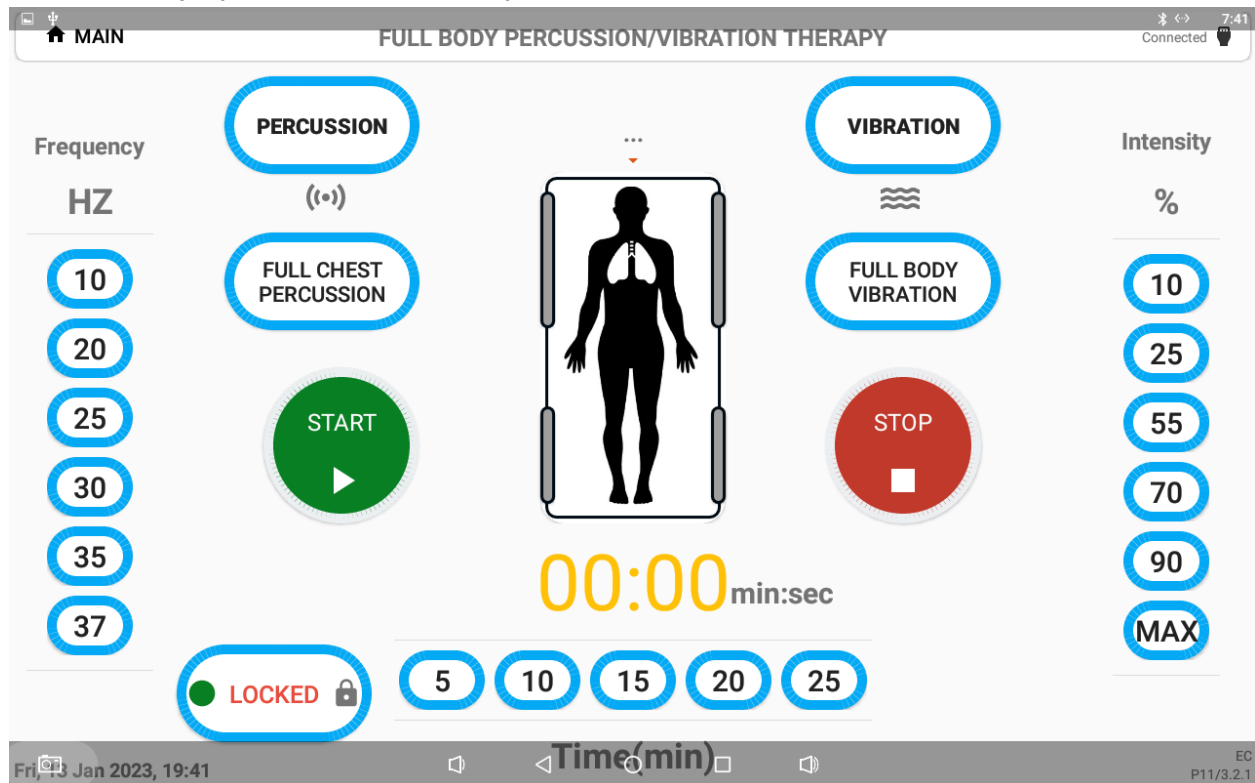


Figure 17-Therapy locked.

To enable and disable press and hold “Locked button” for 4 seconds until the green dot appear or disappear according to the case. When this option is enabled, the system will work until therapy is complete, after that the user can not change any parameter’s until disable locked.

Therapy complete.

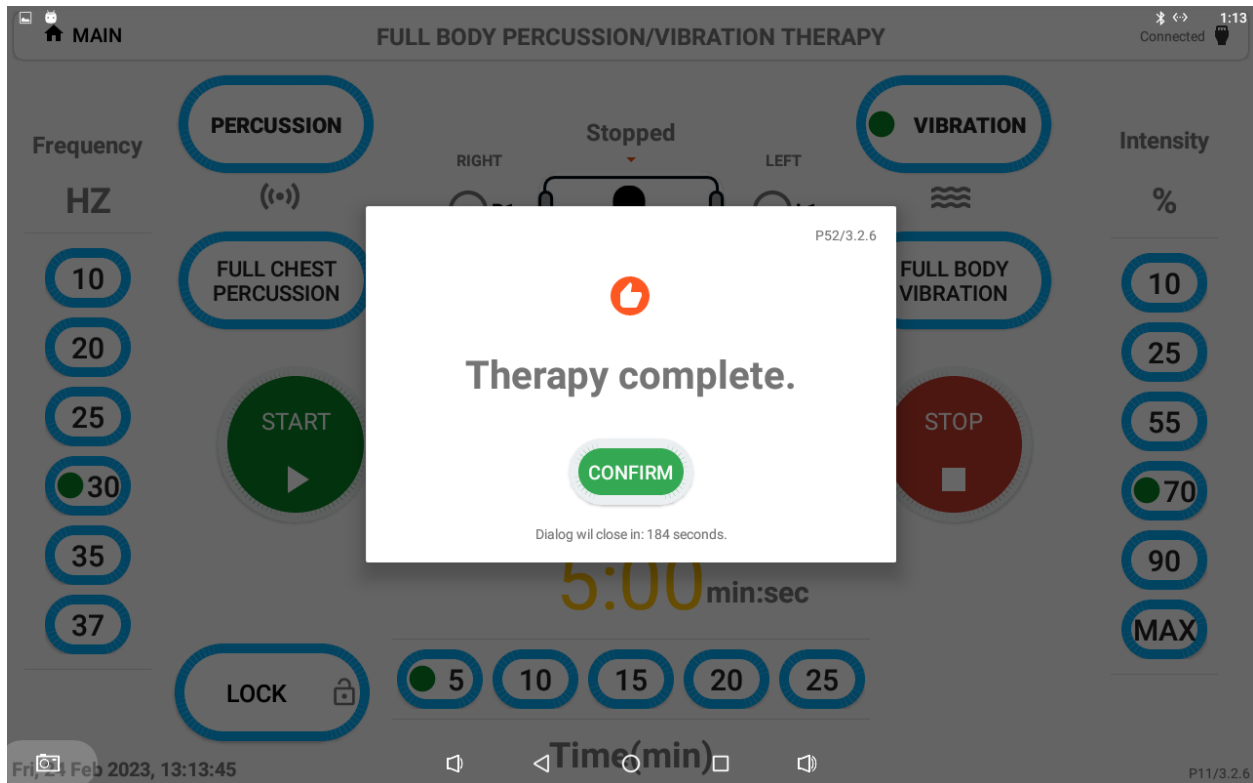


Figure 18-Therapy complete.